

Elekta's purpose is to invent and develop effective solutions for the treatment of cancer and brain disorders. Our goal is to help our customers deliver the best care for every patient. Our oncology and neurosurgery tools and treatment planning systems are used in more than 6,000 hospitals worldwide. They help treat over 100,000 patients every day. The company was founded in 1974 by Professor Lars Leksell, a physician. Today, with its headquarters in Stockholm, Sweden, Elekta employs around 4,000 people in more than 30 offices across 24 countries. The company is listed on NASDAQ OMX Stockholm.

# Monaco® 5.11

## Comprehensive treatment planning

### Personalized Services – Beyond Standard Support


Elekta employs the largest full-time staff of dedicated radiation treatment planning professionals in the industry in the fields of research & development and customer support. Our physics services for beam data modeling will enhance the Monaco user experience and create efficiencies in moving to clinical use status. When you select our solutions, you gain access to our entire team:

- PhD and MS physicists
- PhD mathematicians
- CMDs and RTTs
- Professional software and hardware engineers

#### Support activities include:

- Global live and Web-Support
- Customer training conducted from RT professionals
- Advanced user training

For more information, contact your local Elekta representative.

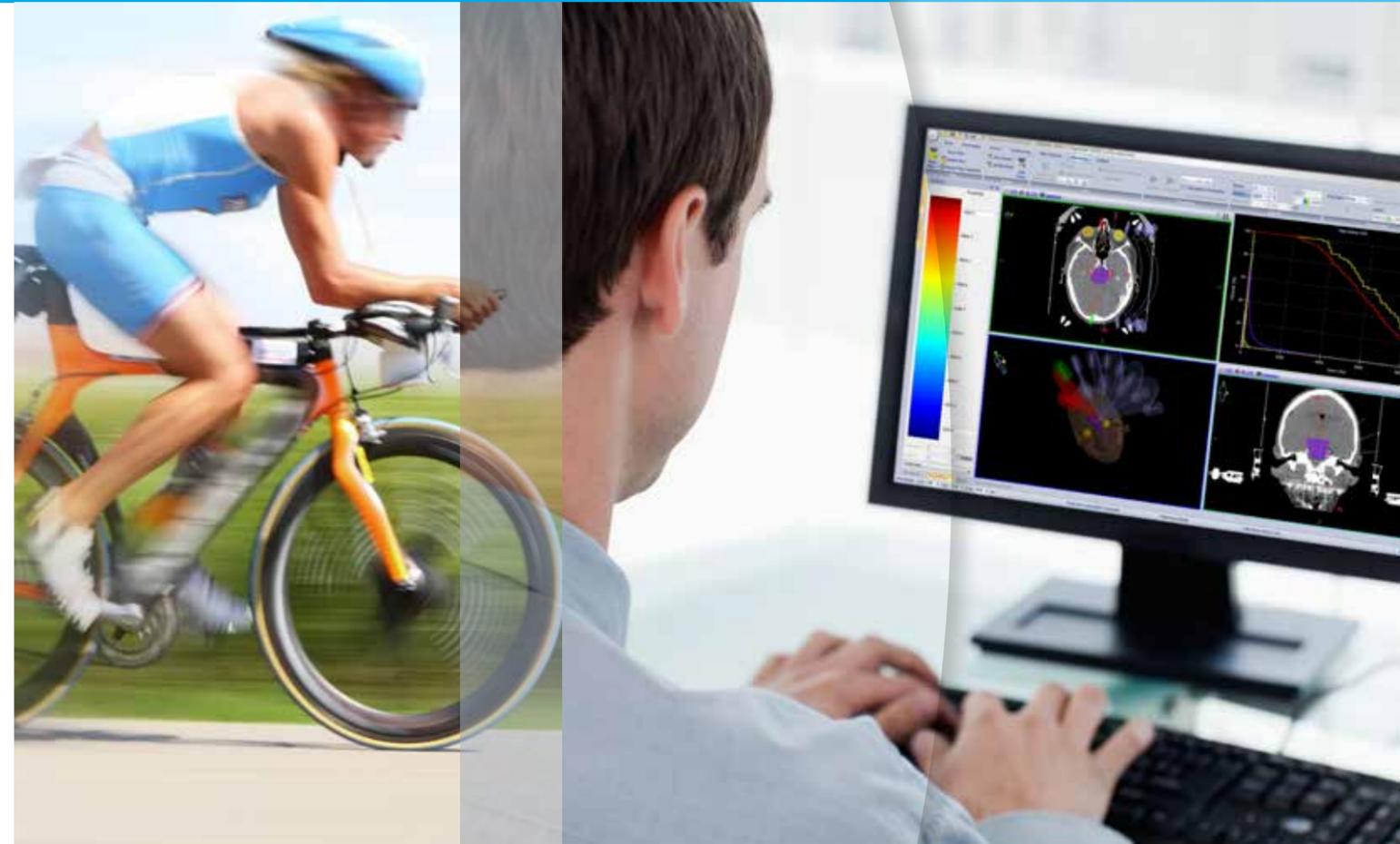

*Accuracy accelerated*



[www.elekta.com](http://www.elekta.com)

*Human Care Makes the Future Possible*

#### Corporate Head Office:

Elekta AB (publ)  
Box 7593, SE-103 93 Stockholm, Sweden  
Tel +46 8 587 254 00  
Fax +46 8 587 255 00  
[info@elekta.com](mailto:info@elekta.com)

#### Regional Sales, Marketing and Service:

##### North America

Tel +1 770 300 9725  
Fax +1 770 448 6338  
[info.america@elekta.com](mailto:info.america@elekta.com)

##### Europe, Middle East, Africa, Eastern Europe, Latin America

Tel +46 8 587 254 00  
Fax +46 8 587 255 00  
[info.europe@elekta.com](mailto:info.europe@elekta.com)

##### Asia Pacific

Tel +852 2891 2208  
Fax +852 2575 7133  
[info.asia@elekta.com](mailto:info.asia@elekta.com)



Art. nr. LPQI0N090520 v6.0 © 2016 Elekta AB (publ). All mentioned trademarks and registered trademarks are the property of the Elekta Group. All rights reserved. No part of this document may be reproduced in any form without written permission from the copyright holder.



# A Complete Planning System with Speed and Accuracy

Monaco, an outstanding planning solution for years, now offers clinics the opportunity to do more, do it faster and do it more accurately than ever before. Monaco now consolidates support for all the most clinically useful delivery techniques in a single system, including flattening filter free and features refinements that bring planning speed and precision to new levels.

## Your Single Source for Planning

With Monaco, you have a comprehensive system to support all major modalities of treatment. It supports 3D conformal radiation therapy, IMRT, VMAT, stereotactic MLC and cones. In addition, it supports Siemens mARC (rIMRT) as well as many of your 4D needs. The simulation package encompasses 4D contouring and the convenient EZ Sketch™ functionality.

## A Higher Standard for Accuracy

Monaco has used the Monte Carlo algorithm – the most accurate dose calculation available – longer than any treatment planning system. Elekta's nearly 40 years of experience in treatment planning have enabled ongoing refinements in this algorithm to ensure the highest possible standard for planning accuracy.

To improve plan quality Monaco offers Segment Shape Optimization™, which strives to match optimized plan doses while meeting linac deliverability. In addition, the unique biological modeling capabilities of Monaco apply defined clinical objectives, maximizing confidence in plan assessment and selection.

Multi-criterial optimization assures organs at risk are spared to the highest degree possible while maintaining target coverage. Real time interaction during and after optimization provides the flexibility to precisely tune the plan results as needed.

## Satisfying Your Need for Planning Speed

Monaco now provides a number of enhancements that accelerate your workflow. These include a better organized user interface configuration that improves access to tools, and a planning control toolbar that provides immediate access to beam or structure spread sheets, prescription, Dose reference points (DRP) and IMRT constraints.

In addition, compared to previous versions, Monte Carlo control point dose calculations are faster\*. The addition of the Collapsed Cone algorithm with an available Graphics Processing Unit (GPU) card helps plan calculation efficiency.

Monaco now features a sleek and easy-to-follow user interface to guide you to advanced features intended to create accurate plans for you and your patients. Utilize templates to gain greater planning efficiency. Create highly conformal treatment plans quickly and easily, and share with other centers.

## Why Monaco®?

- Comprehensive, single box solution
- Personalized packages to meet your clinical needs with future upgrade opportunities to Elekta Software
- Optimized code utilizing gold-standard Monte Carlo algorithm offers accuracy and speed
- Stereotactic treatment planning support
- Unique biological modeling capabilities
- User interface tools accelerate planning
- Leading support and service team to advise you at every stage and level of treatment planning

\*Data on file

# Bring Monaco into Your Clinical Practice

## Enhance Planning Performance

Choose from a broad suite of planning tools and range of robust dose calculation algorithms. Clinicians can select the tools for the technique to optimize radiation therapy delivery for each patient. With a new, intuitive and user-friendly interface, Monaco allows users to quickly access important plan information including beam or structure spreadsheets, prescription, DRP, and IMRT constraints.

## 3D Conformal

With the addition of 3D conformal treatment planning capabilities, Monaco includes customizable

templates that support efficient plan generation. By using a Graphics Processing Unit (GPU) for the collapsed-cone algorithm, Monaco increases the speed and performance of 3D dose calculations even further. And, Monaco supports a range of modalities to cater to the clinic's needs for expanding the treatment planning capabilities, including support of wedges, bolus and VMC++ Electron Monte Carlo. Use Forward Planning to achieve better dose uniformity than can be obtained with a single-port 3D beam.

## IMRT

Monaco integrates innovative

biological cost functions with multicriterial constrained optimization, a powerful leaf-sequence optimizer and a robust Monte Carlo dose calculation algorithm to create the most advanced IMRT planning solution available today. Monaco also can be used effectively in most clinical environments. As a vendor neutral IMRT planning platform, Monaco supports all major linear accelerators and connects to any record-and-verify information system, including Elekta's MOSAIQ® Oncology Information System. Monaco includes fast pencil beam dose engine for "ideal" fluence optimization.




## VMAT

Monaco's Volumetric Modulated Arc Therapy (VMAT) functionality can optimize single or multiple non-coplanar arcs simultaneously, providing the flexibility and control needed for more complex treatment plans. Arc plans can be delivered with a single button push at the linear accelerator console. Gantry directions are automatically sorted and all control points are seamlessly integrated into a single deliverable arc sequence. Monaco offers the XVMC Monte Carlo dose engine, for electron and photon, for a continuous arc calculation as a single beam, rather than dose

approximations that occur with many discrete (control point) gantry angle positions. Monaco also offers the Dynamic Conformal Arc feature.




## Integration

Monaco connects to Elekta's MOSAIQ for a truly optimized clinical workflow and also supports all major OIS and Linac connectivity. Monaco fully integrates Elekta's superior contouring and fusion capabilities, including the ability to contour 4D CT data, with the option of providing ABAS on the same hardware for a seamless transition from delineation to planning.

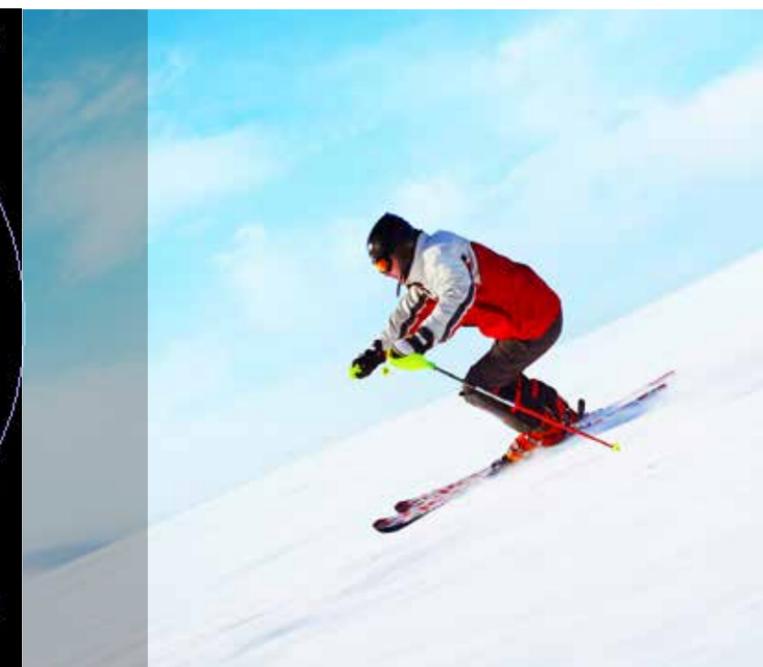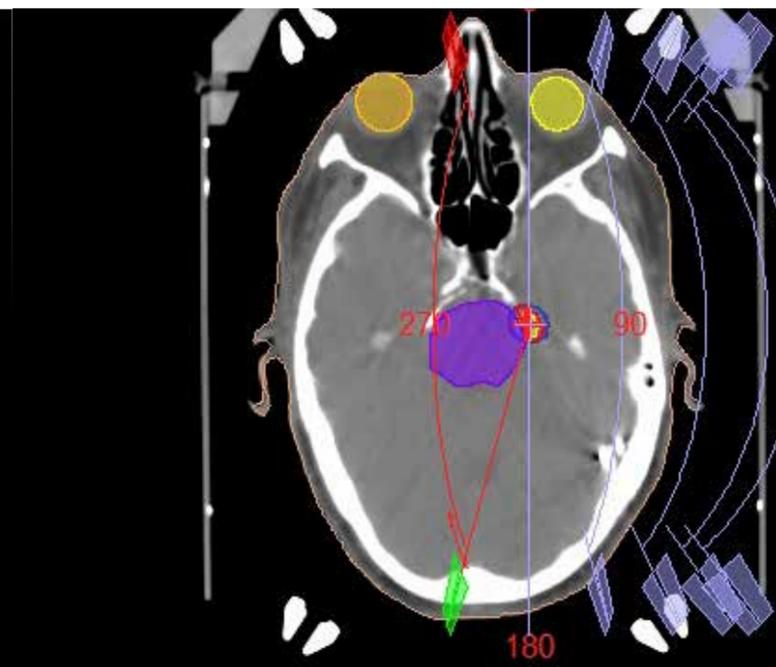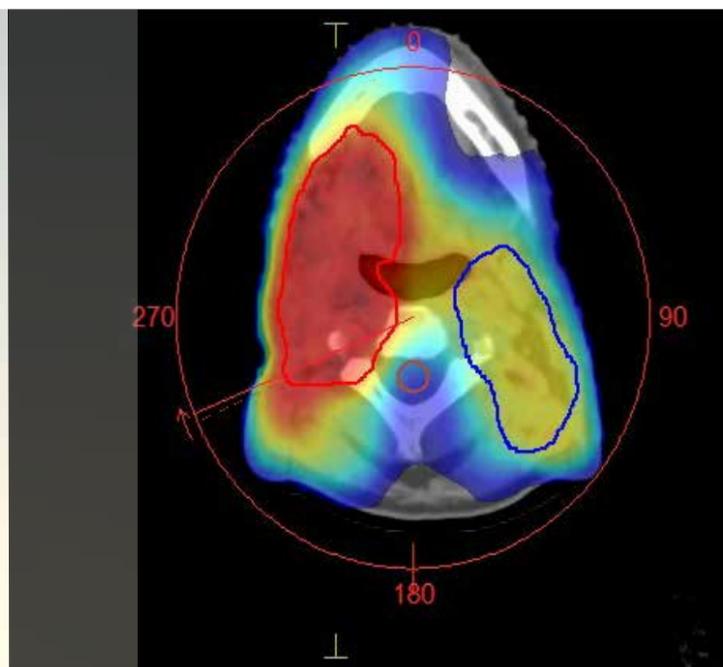

# Monaco Advances

Monaco® features advanced concepts to enable complete planning for all techniques; complete stereotaxy planning for frame-based, frameless, APEX™ micro-leaf MLC, and cones; contouring for 3D and 4D data and more. These unique features drive much of the functionality for Monaco.

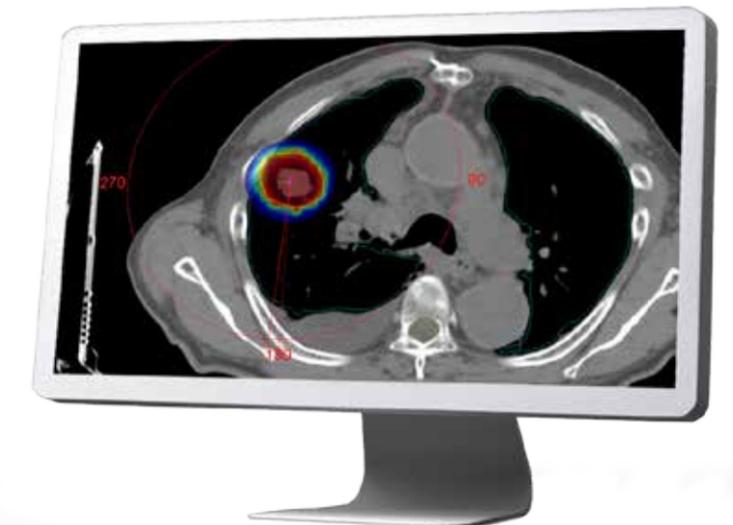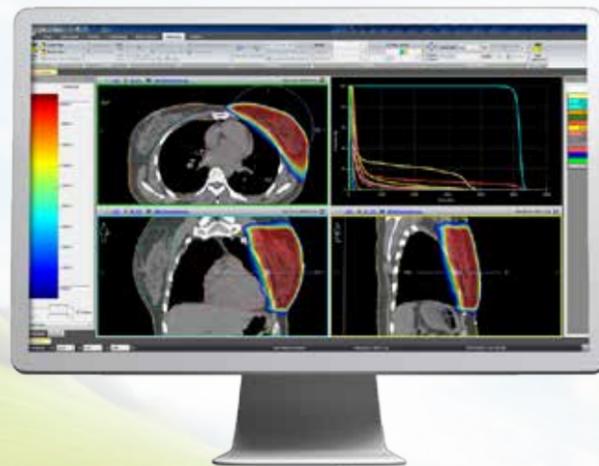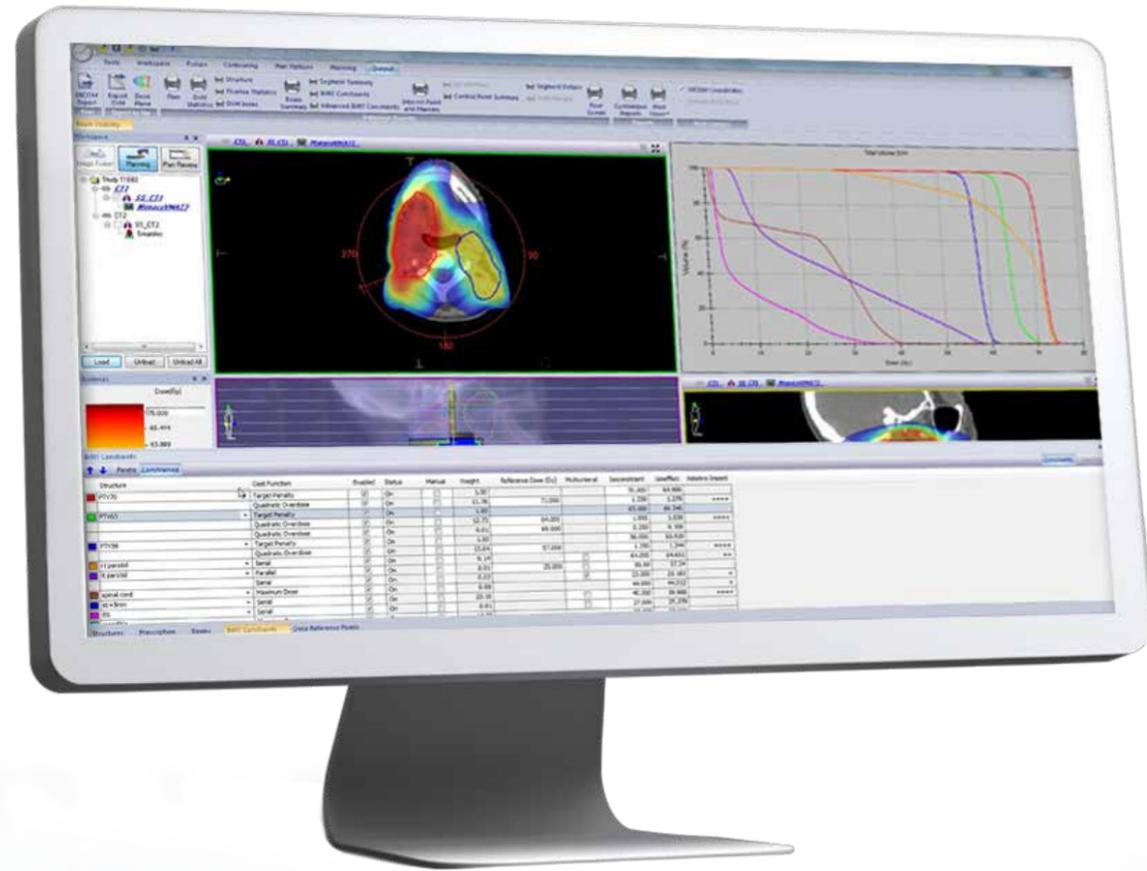

## Biological Modeling

- Wide range of optimization cost functions, including dose-volume-based and biological cost functions
- Equivalent Uniform Dose and iso-effective volumes account for tissue response to dose per fraction and volume effect of organs
- Complete control of the DVHs for healthy tissue and tumor volumes

## Multi-criterial Goals

- Multi-criterial optimization offers precise adjustments generated in real-time during optimization
- Automated tightening of constraints to optimize normal tissue sparing
- Fine tune the plan during the process and after the final calculation with constraint weight edits and a fast update

## Sensitivity Analysis

- Unique sensitivity analysis tool easily resolves conflicts between target dose objective and dose limiting constraints

Optimization evaluation tool shows interdependence of each constraint and objective on target dose coverage

## Monte Carlo Dose Engine

- XVMC Monte Carlo dose engine – for accurate dose calculation
- Fast pencil beam dose engine for “ideal” fluence optimization – specifically created for IMRT
- Includes both electron and photon beams

## Constrained Optimization

- Ensures that OAR dose limits are satisfied
- Changes to any constraint in the prescription don't affect the other normal tissue constraints, only the dose to target(s)
- Enhanced voxel controls guide dose gradients between structures, making further “technical” volumes unnecessary

## Sequencing

### Segment Shape Optimization™

- Proprietary approach to smoothing and clustering segments, then optimizing beam weights and shapes enables clinicians to improve dose conformity, plan quality and delivery efficiency
- Reduce the number of segments, shortening treatment times
- Powered by Monte Carlo, develop plans that improve sparing of OARs and dose conformity

### Stereotactic Planning

- Combine Monaco with Versa HD® for the system accuracy you need for the best, high-definition stereotactic plans and efficient deliveries.
- Dynamic Conformal Arc (DCA) Therapy enables clinicians to deliver highly conformal stereotactic plans using the MLC to dynamically conform around a target as the beam rotates around the patient
- Create multiple arcs in a non-coplanar fashion to render more conformal plans
- Used in conjunction with micro-MLCs, such as Elekta's Apex™ MLC, it can yield very conformal plans using DCA or VMAT
- Supports stereotactic cone treatments for static and arc fields

### MR Planning

- Create radiotherapy treatment plans on transverse MR images using forced electron densities for the dose calculation