



MOSAIQ IQ Script: Create Assessment by Appointment Activity

Create and assign preference to present a COVID-19 assessment when an appointment is scheduled, ended and/or queued


- 1) Click **File > System Utilities > IQ Engine Configuration**
- 2) Click **Scripts > Operations and Events > General Charting > Appointment Script**
- 3) Click **Create Assessment by Appointment Activity**
- 4) Under preferences, click **New**
- 5) Select **Assessment/Tab View**

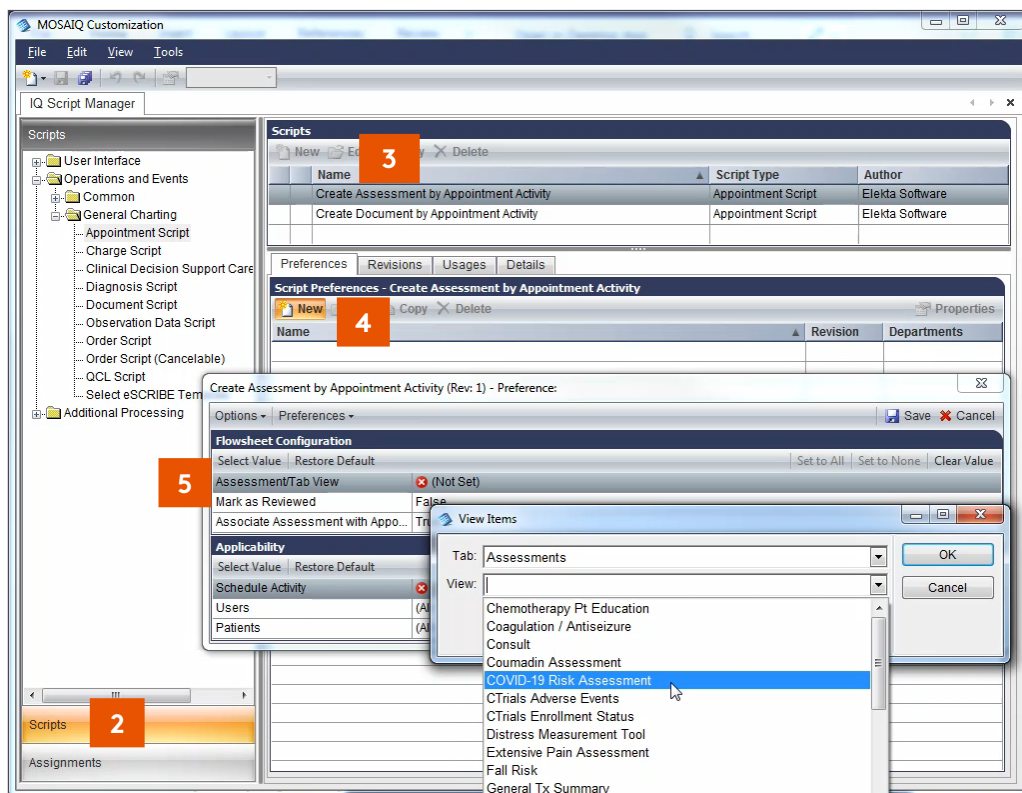

- 7) Select one or more **Schedule Activities**
- 8) Optionally, set **Users** or **Patients** (for testing purposes only)

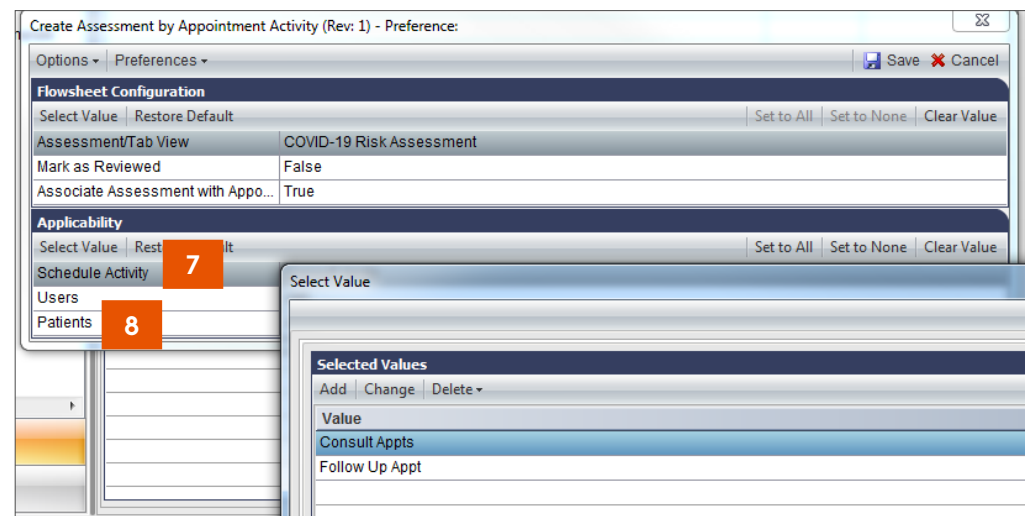

- 9) Name preference, select all relevant departments and click **OK** to save

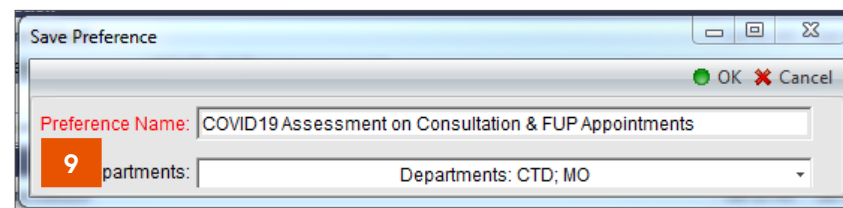


- 6) Optionally, set **Mark as Reviewed** to true and/or **Associate Assessment with appointment** to False


MOSAIQ IQ Script: Create Assessment by Appointment Activity continued

Create and assign preference to present a COVID-19 assessment when an appointment is scheduled, ended and/or queued

- 10) Click **Assignments > Operations and Events > General Charting > Appointment Script > Appointments**
- 11) Click one of **Appointment Ended, Appointment Scheduled** or **Patient Queued**
- 12) Under Script Assignments, click **New**
- 13) Select **Create Assessment by Appointment Activity** script

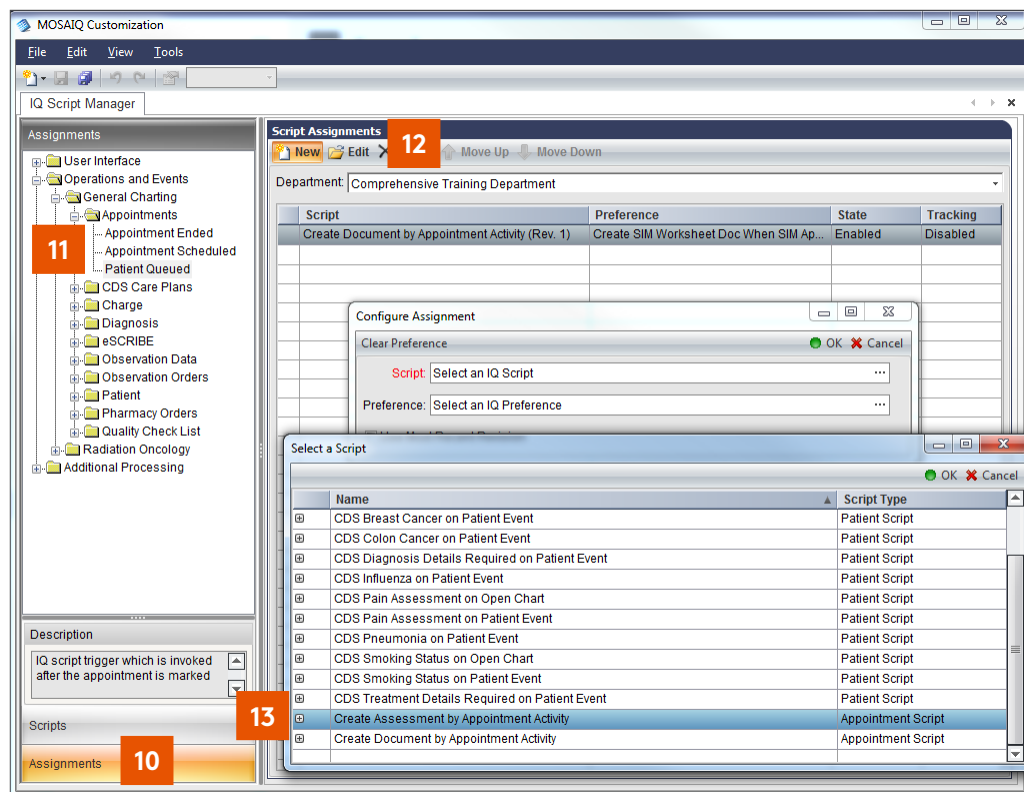

- 14) Select **Preference**

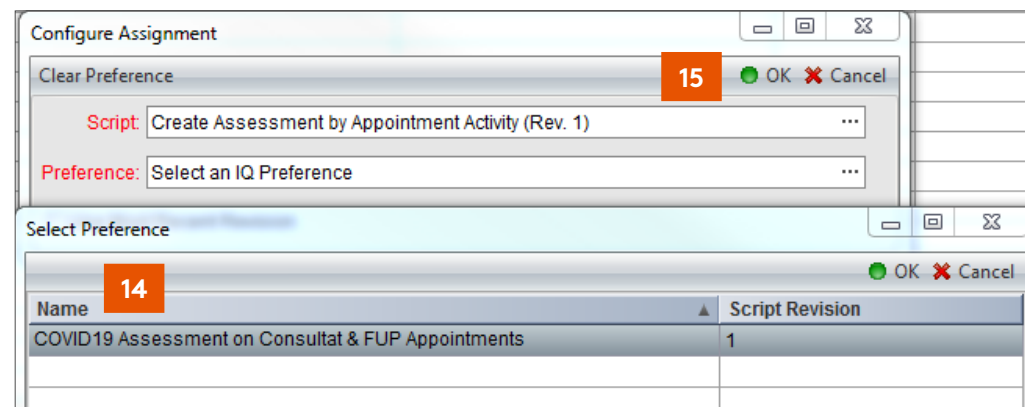

15. Click **OK** to save configuration assignment
16. Remember to repeat for all relevant departments

