
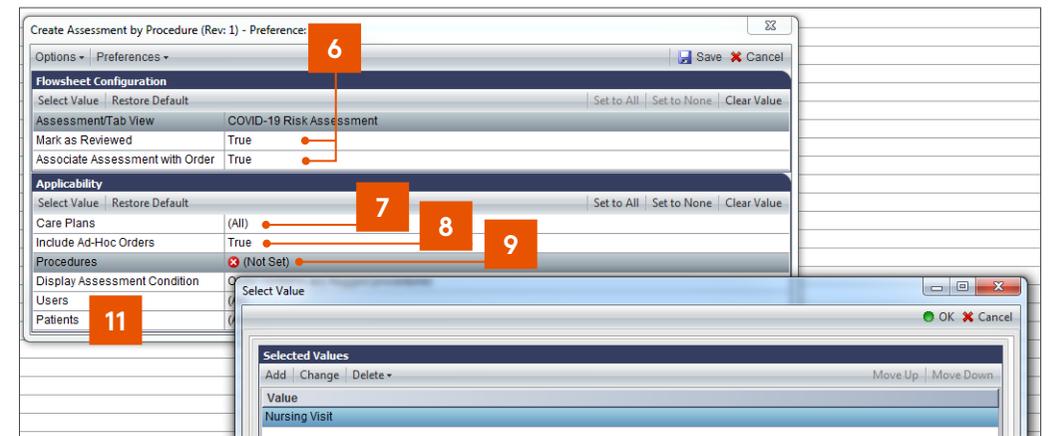


# MOSAIQ IQ Script: Create Assessment by Procedure


Create and assign preference to present a COVID-19 assessment when an observation order (procedure) is being approved or generated approved through Quick RX

- 1) Click **File > System Utilities > IQ Engine Configuration**
- 2) Click **Scripts > Operations and Events > General Charting > Order Script**
- 3) Click **Create Assessment by Procedure**
- 4) Under preferences, click **New**
- 5) Select **Assessment/Tab View**

- 6) Optionally, set **Mark as Reviewed** to true and **Associate Assessment with Order** to false
- 7) Optionally select one or more **Care Plans**
- 8) Optionally set **Include Ad-Hoc Orders** to false
- 9) Select one or more **Procedures**
- 10) If selecting more than one procedure, optional update **Display Assessment Condition** to **Order contains any flagged procedure**
- 11) Optionally, set **Users** or **Patients** (for testing purposes only)

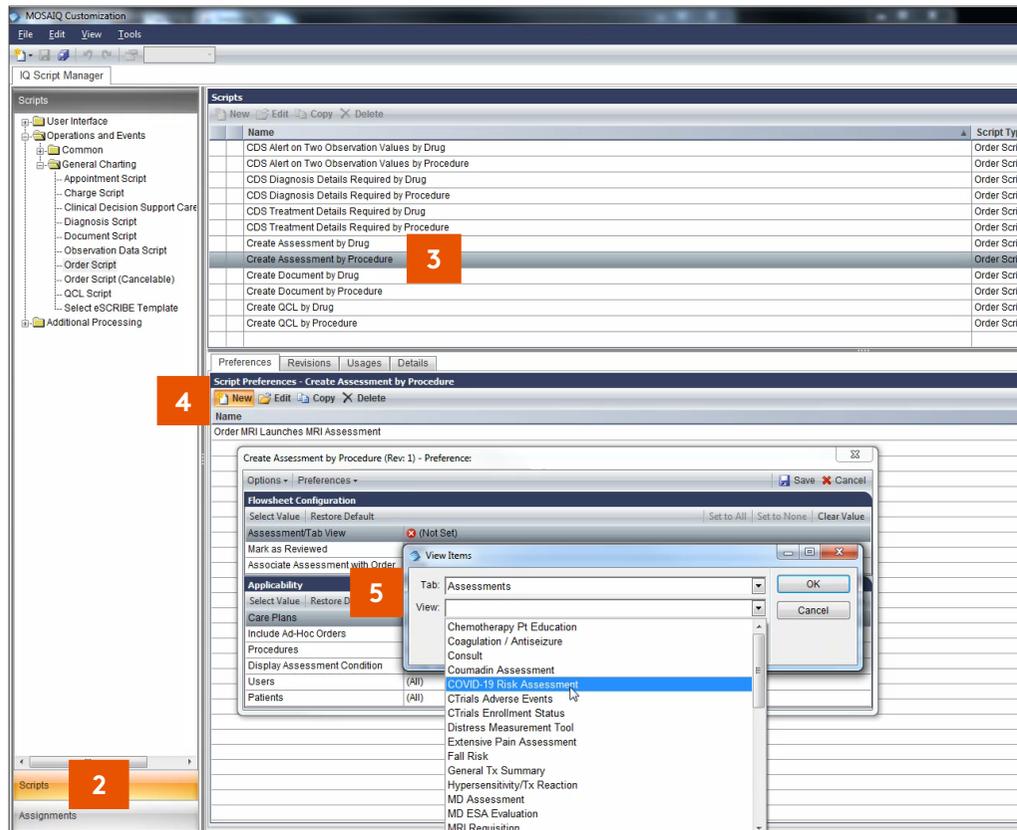


- 12) Name preference, select all relevant departments and click **OK** to save

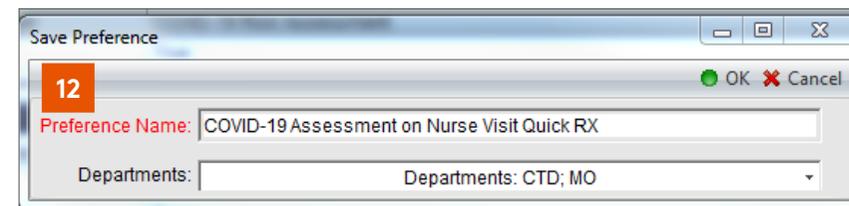


# MOSAIQ IQ Script: Create Assessment by Procedure continued

Create and assign preference to present a COVID-19 assessment when an observation order (procedure) is being approved or generated approved through Quick RX

- 13) Click **Assignments > Operations and Events > General Charting > Observation Orders**
- 14) Click **Order Approved**
- 15) Under Script Assignments, click **New**
- 16) Select **Create Assessment by Procedure** script

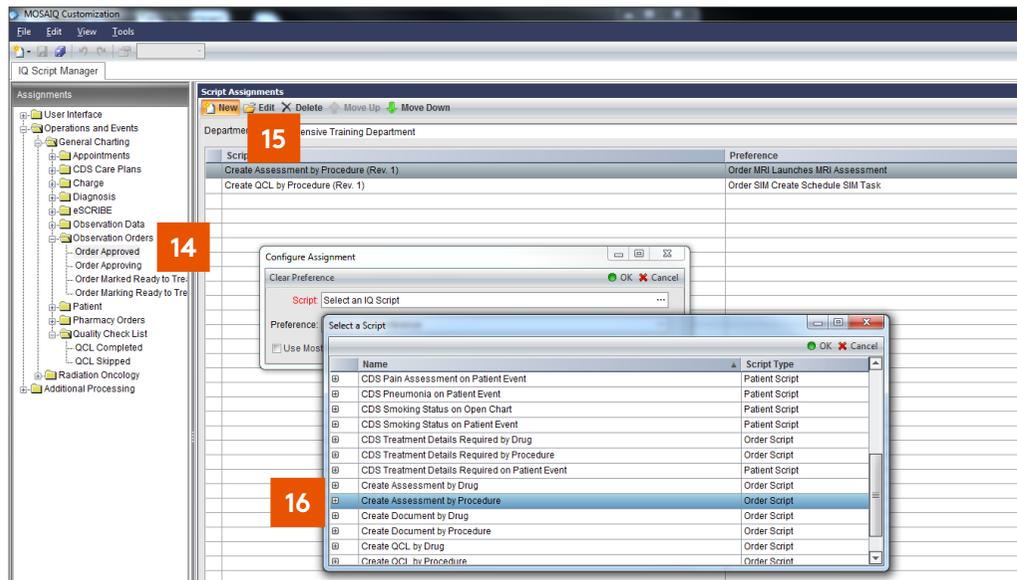

- 17) Select **Preference**

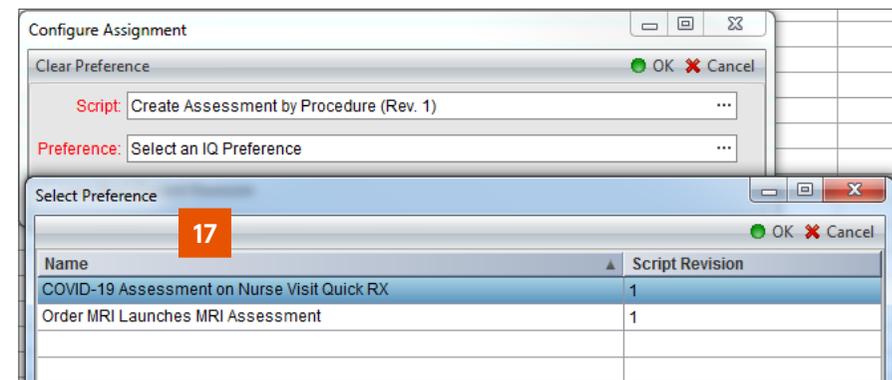

- 18) Click **OK** to save configuration assignment
- 19) Remember to repeat for all relevant departments

