



# MOSAIQ IQ Script: Create Assessment by QCL Task

Create and assign preference to present a COVID-19 assessment when a QCL Task is completed or skipped


- 1) Click **File > System Utilities > IQ Engine Configuration**
- 2) Click **Scripts > Operations and Events > General Charting > QCL Script**
- 3) Click **Create Assessment by QCL Task**
- 4) Under preferences, click **New**
- 5) Select **Assessment/Tab View**

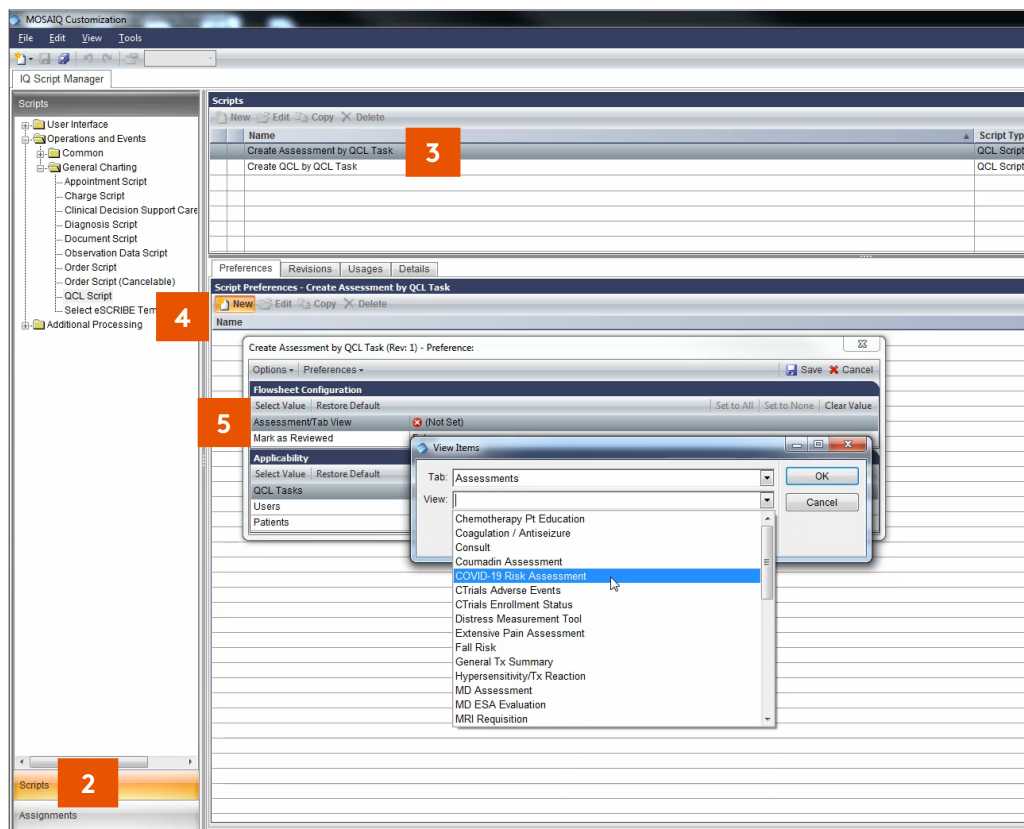

- 6) Optionally, set **Mark as Reviewed** to true
- 7) Select one or more **QCL Tasks**
- 8) Optionally, set **Users** or **Patients** (for testing purposes only)

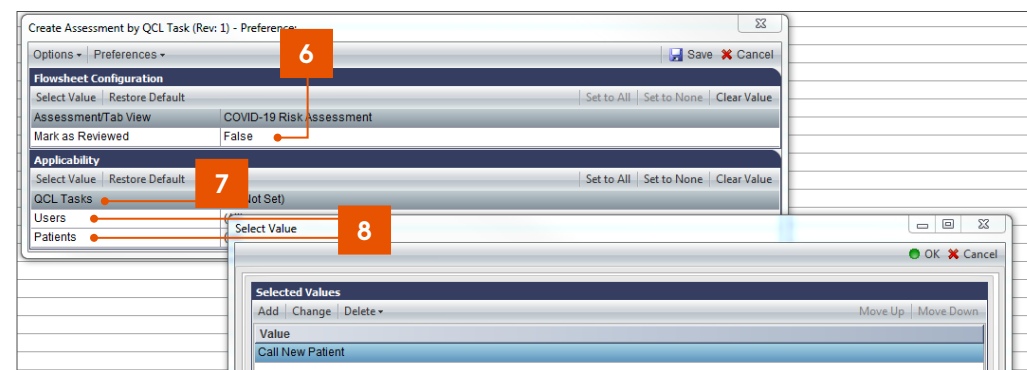


- 9) Name preference, select all relevant departments and click **OK** to save

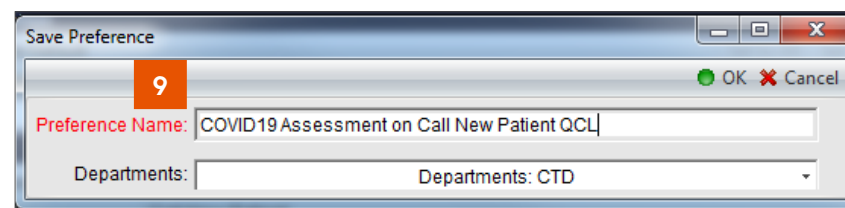


# MOSAIQ IQ Script: Create Assessment by QCL Task continued

Create and assign preference to present a COVID-19 assessment when a QCL Task is completed or skipped

- 10) Click **Assignments > Operations and Events > General Charting > Quality Check List**
- 11) Click one of **QCL Completed** or **QCL Skipped** (most commonly completed)
- 12) Under Script Assignments, click **New**
- 13) Select **Create Assessment by by QCL Task** script

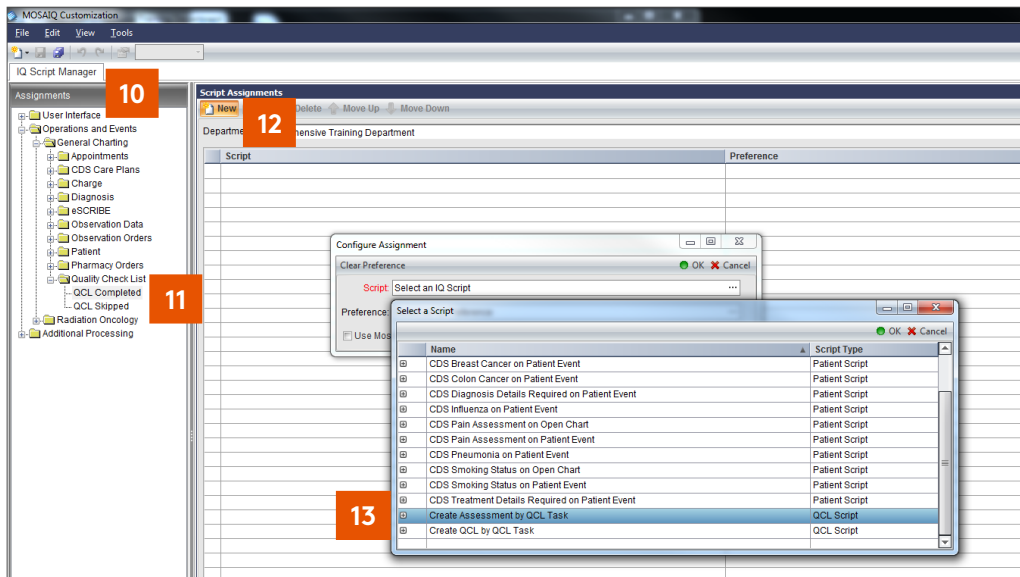

- 14) Select **Preference**

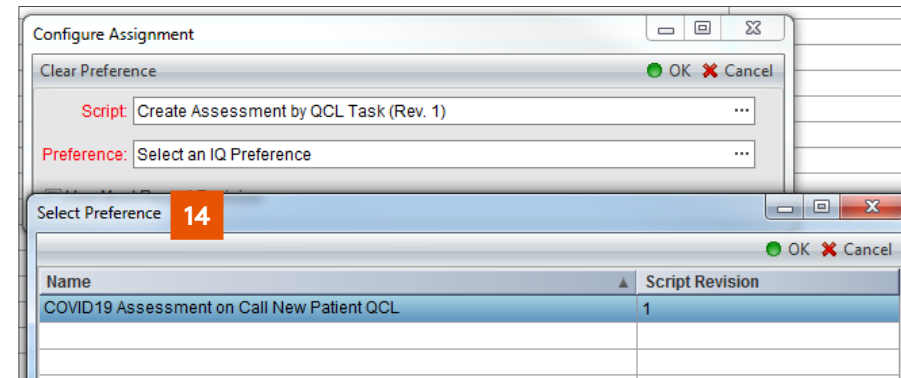

- 15) Click **OK** to save configuration assignment
- 16) Remember to repeat for all relevant departments

