
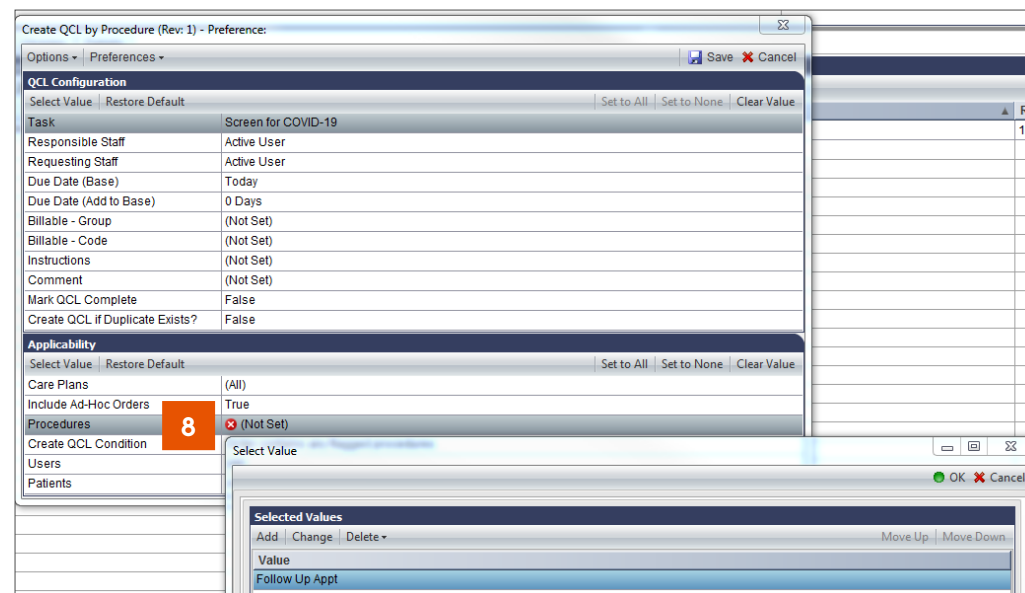


# MOSAIQ IQ Script: Create QCL by Procedure

Create and assign preference to create a QCL task when an observation order (procedure) is being approved or generated approved through Quick Rx

- 1) Click **File > System Utilities > IQ Engine Configuration**
- 2) Click **Scripts > Operations and Events > General Charting > Order Script**
- 3) Click **Create QCL by Procedure**
- 4) Under preferences, click **New**
- 5) Select **Task**


- 6) Optionally, set **Responsible Staff, Requesting Staff, Due Date (Base), Due Date (Add to Base), Billable Group, Billable Code, Instructions, Comment, Mark QCL Complete, Create QCL if Duplicate Exists?**
- 7) Optionally select one or more **Care Plans** and/or **Include Ad-Hoc Orders** to false
- 8) Select one or more **Procedures**
- 9) If selecting more than one procedure, optionally update **Create QCL Condition** to **Order contains any flagged procedure**

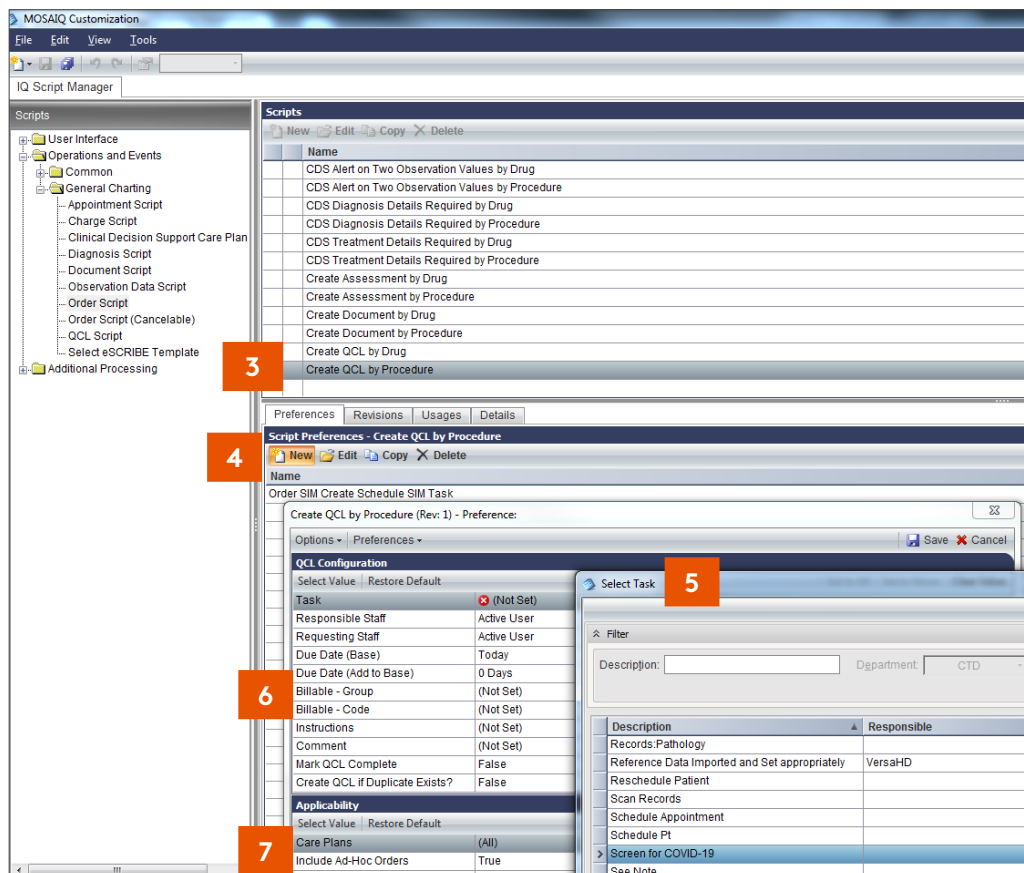


# MOSAIQ IQ Script: Create QCL by Procedure continued

Create and assign preference to create a QCL task when an observation order (procedure) is being approved or generated approved through Quick Rx


- 10) Optionally, set **Users** or **Patients** (for testing purposes only)
- 11) Name preference, select all relevant departments and click **OK** to save

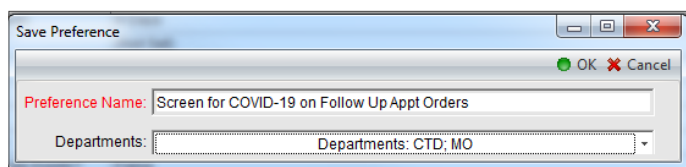

- 12) Click **Assignments > Operations and Events > General Charting > Observation Orders**
- 13) Click **Order Approved**
- 14) Under Script Assignments, click **New**
- 15) Select **Create QCL by Procedure** script

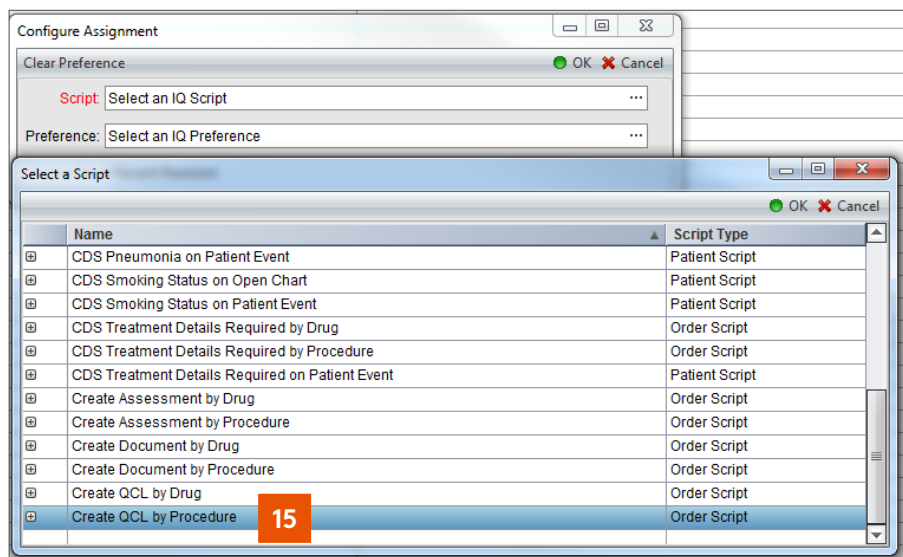

- 16) Select **Preference**

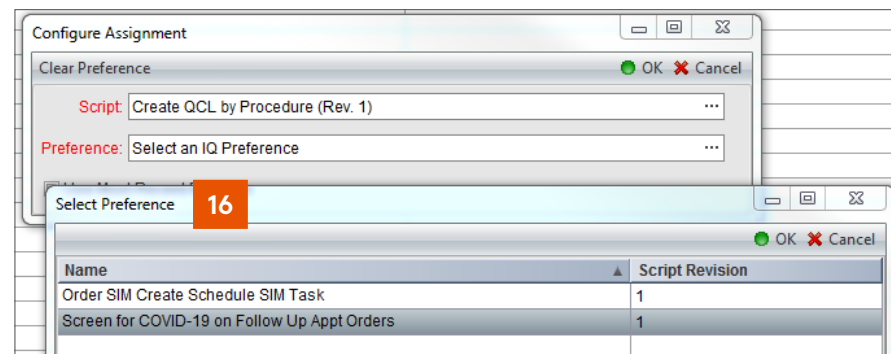

- 17) Click **OK** to save configuration assignment
- 18) Remember to repeat for all relevant departments

Script	Preference	State	Tracking
Create Assessment by Procedure (Rev. 1)	Order MRI Launches MRI Assessment	Enabled	Disabled
Create QCL by Procedure (Rev. 1)	Order SIM Create Schedule SIM Task	Enabled	Disabled
Create Assessment by Procedure (Rev. 1)	COVID-19 Assessment on Nurse Visit Quick RX	Enabled	Disabled
Create QCL by Procedure (Rev. 1)	Screen for COVID-19 on Follow Up Appt Orders	Enabled	Disabled