


# MOSAIQ IQ Script: Create QCL on Observation Value


Create & assign preference to create a QCL task when a response on a COVID-19 assessment (observation value) requires further action

- 1) Click **File > System Utilities > IQ Engine Configuration**
- 2) Click **Scripts > Operations and Events > General Charting > Observation Data Script**
- 3) Click **Create QCL on Observation Value**
- 4) Under preferences, click **New**
- 5) Select **Task**

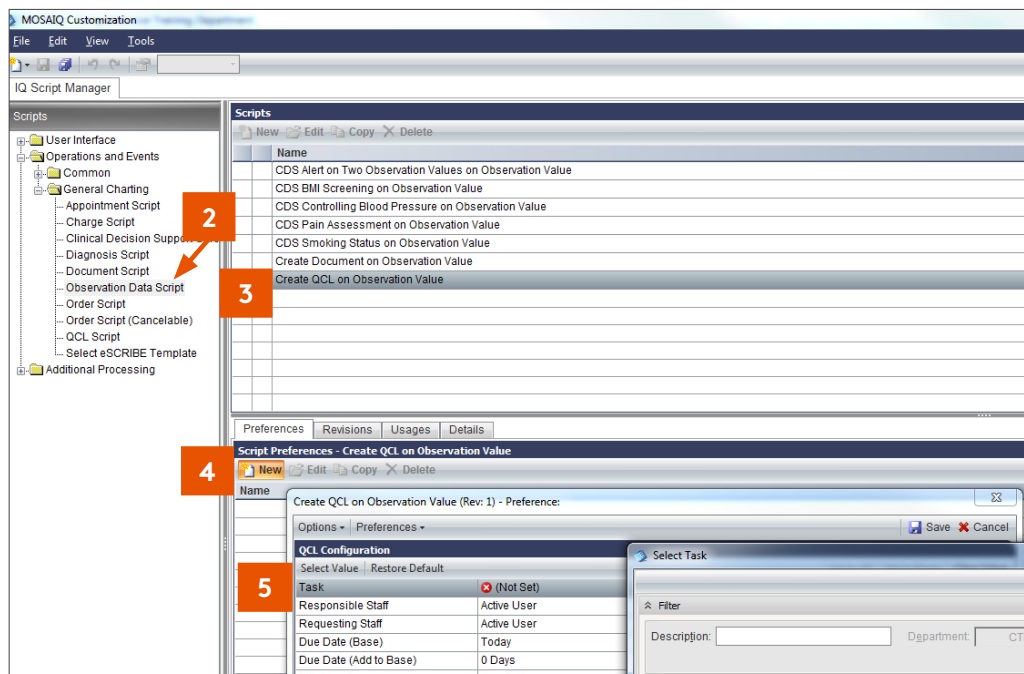


- 6) Optionally, set **Responsible Staff, Requesting Staff, Due Date (Base), Due Date (Add to Base), Billable Group, Billable Code, Instructions, Comment, Mark QCL Complete, Create QCL if Duplicate Exists?**

- 7) Optionally, set **Which Diagnosis, Diagnosis Code, View**
- 8) Select **Data Items with Thresholds**

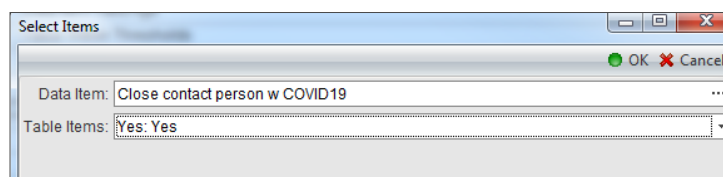


- 9) Optionally set **Valid Age, Data Item Filter, Data Item Filter Mode**
- 10) Optionally, set **Users or Patients** (for testing purposes only)
- 11) Name preference, select all relevant departments and click **OK** to save

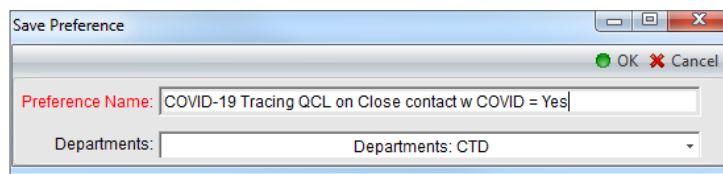


## MOSAIQ IQ Script: Create QCL on Observation Value continued

Create & assign preference to create a QCL task when a response on a COVID-19 assessment (observation value) requires further action

- 12) Click **Assignments > Operations and Events > General Charting > Observation Data**
- 13) Click one of **Observation Data Created** or **Observation Data Reviewed**
- 14) Under Script Assignments, click **New**
- 15) Select **Create QCL on Observation Value** script

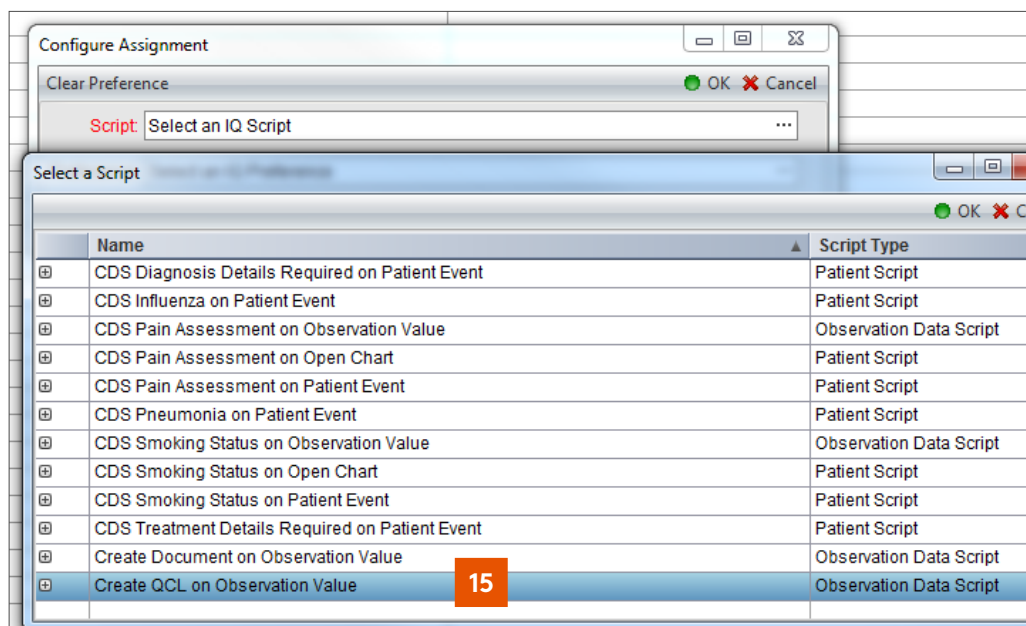

- 16) Select **Preference**

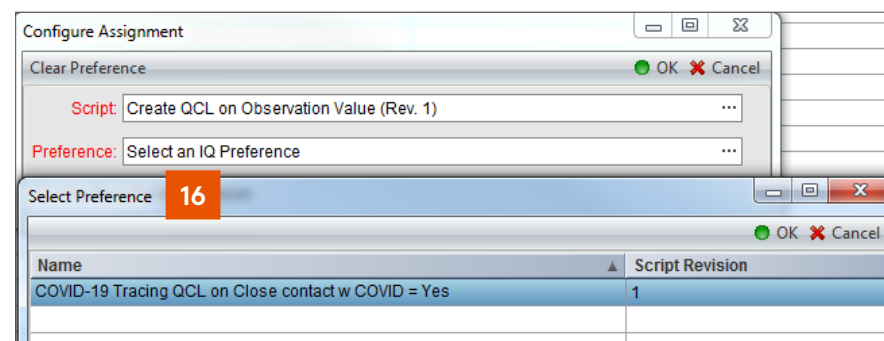

- 17) Click **OK** to save configuration assignment
- 18) Remember to repeat for all relevant departments

